

SAFE PERMANENT MAKE-UP

bio**E**volution





BioEvolution, in cooperation with PowerLiner.pl, based on over 20 years of experience, has launched the production of high-class permanent makeup professional devices, hand pieces and cartridges in Poland. Thanks to our own production line we can keep control over quality and precision at every stage of product development. Modern and ergonomic design, state-of-the-art materials, extraordinary attention to detail and innovative solutions are the characteristic features of the new concept.

## TECHNOLOGY

Power  Liner  
TECHNOLOGY INSIDE


# ACUPUNCTURE IN THE CARTRIDGE

secure  
CARTRIDGE


SECURE cartridges are the first in the world that are available in two varieties: with acupuncture needles and traditional needles. The innovative use of acupuncture needles in cartridges makes using the acupuncture technique not only extremely precise, but also safe.

## POWER AND PRECISION

BioEvolution devices are characterized by the high needle power, and precision - regardless of the density and hardness of the skin. This minimises bleeding and down-time for the client. Thanks to the precise nature of the needle, the pigment line is perfectly consistent, and not jagged. The innovative POWER LINER technology used in the hand piece ensures continuous frequency, power, and predictability - regardless of changing skin structures. The benefits of this feature are more pigment in the skin, at a more consistent rate - meaning a faster treatment, quicker healing time, and a more precise treatment.

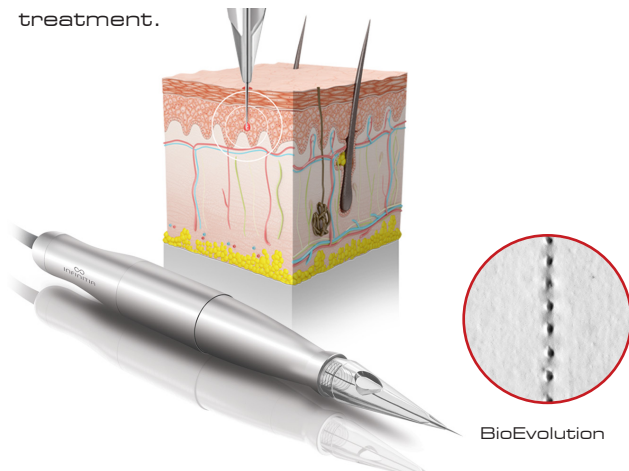


Traditional permanent makeup devices press down on the skin under the pressure of the needle during treatment. This can cause the skin to 'pull up' when the needle is moved away from the skin. Insufficient power and frequency of the needle also causes problems for the needle overcoming the natural resistance of the surface of the skin. This issue has been overcome with the new system, and we no longer see issues with surface skin damage, or heavy bleeding in the treatment area - which also causes pigment flush from the skin.

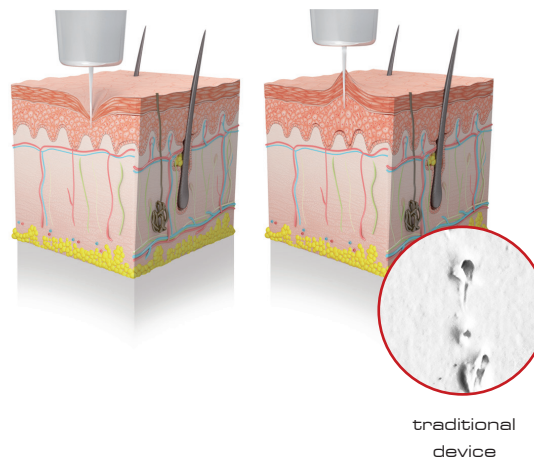

## MORE COLOUR

The specially designed shape of the SECURE needle tip enables the insertion of up to 50% more pigment during the pass. Needles are manufactured with a 1 micron tolerance for unmatched dexterity and precision. This allows the technician to carry out a more detailed treatment, with more pigment, less flushing, and less downtime for the client

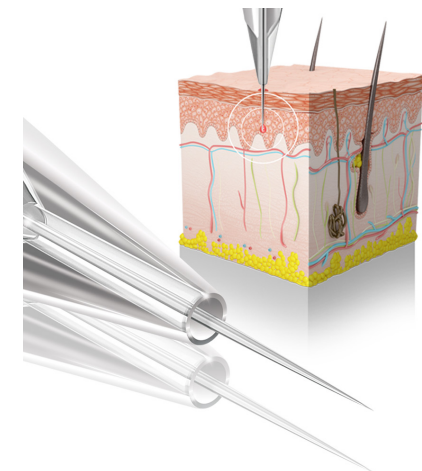




# HANDPIECE

  
INFINITIA

The INFINITIA permanent makeup handpiece has been designed for the most demanding micropigmentation specialist. The innovative drive system ensures the highest precision of the needle. Thanks to the modern and ergonomic design - even long work with the hand piece is not tiring. The highest quality anodized aluminium was used for the production of the device. Thanks to the properly selected weight and precision of workmanship the hand piece operation is very quiet.

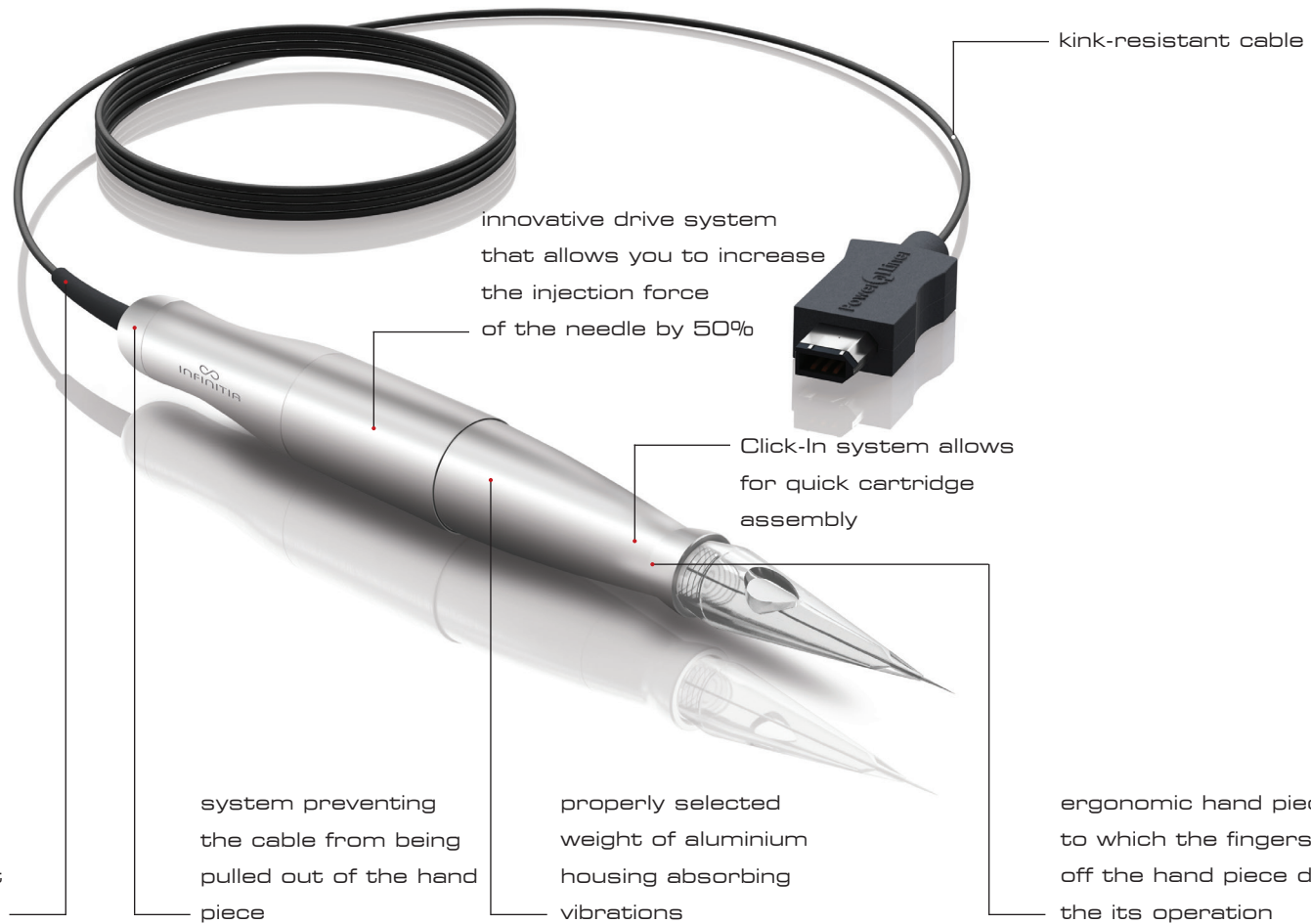

the strain relief protects the most critical section of the cable against kinking

system preventing the cable from being pulled out of the hand piece

innovative drive system that allows you to increase the injection force of the needle by 50%

Click-In system allows for quick cartridge assembly

properly selected weight of aluminium housing absorbing vibrations

ergonomic hand piece, thanks to which the fingers do not slide off the hand piece during the its operation

kink-resistant cable

**Power Liner**  
TECHNOLOGY INSIDE



# VENA

- Needle operational frequency: from 50 to 150 injections per second (adjustment every 10).
- The ergonomic aluminium hand piece is characterized by very quiet operation without any oscillation or vibration. It does not overheat even during a prolonged treatment.
- The SECURE disposable needle cartridge system / the PROTECT system (patent pending), minimizes the possibility of infection without the need to sterilize individual parts in an autoclave.
- The Click-In system allows for quick cartridge assembly.
- Smooth needle length adjustment.
- Control panel: glass and aluminium.
- Foot switch.
- Leather, hand-sewn bag.

#### Device dimensions:

length	148 mm
width	148 mm
height	30 mm

#### Weight:

Device	608 g
hand piece stand	148 g
hand piece with a cable	106 g
hand piece without a cable	78 g
foot switch	112 g

#### Device parameters:

Power supply	230 AC/ 24 DC
Power	10 W
Max. power consumption	1.0 A

#### Hand piece motor parameters:

Power supply	12 V DC
Power	2.57 W
Torque	3.2 Nm
Max. speed	12 000 rpm





PowerLiner  
TECHNOLOGY INSIDE



made in Poland

[www.bioevolution.com](http://www.bioevolution.com)



LABO

- Needle operational frequency: 60, 90, 120 injections per second.
- The ergonomic aluminium hand piece is characterized by very quiet operation without any oscillation and vibration. It does not overheat even during prolonged treatment.
- The SECURE disposable needle cartridge system / the PROTECT system (patent pending), minimizes the possibility of infection without the need to sterilize individual parts in an autoclave.
- The Click-In system allows for quick cartridge assembly.
- Smooth needle length adjustment.
- Control panel: glass and aluminium.

**Device dimensions:**

length	148 mm
width	89 mm
height	23 mm

**Weight:**

Device	302 g
hand piece stand	148 g
hand piece with a cable	106 g
hand piece without a cable	78 g

**Device parameters:**

Power supply	230 AC/ 24 DC
Power	10 W
Max. power consumption	1.0 A

**Hand piece motor parameters:**

Power supply	12 V DC
Power	2.57 W
Torque	3.2 Nm
Max. speed	12 000 rpm



Power  Liner  
TECHNOLOGY INSIDE

made in Poland

[www.bioevolution.com](http://www.bioevolution.com)



## NEEDLE CARTRIDGES

Cartridge - an integrated system of a needle with an adaptor ensures the highest standards of safety and hygiene. The needle in the plastic housing is kept sterile clean. The internal filter, while the cartridge is used, minimizes the amount of colour and blood returning to the hand piece. Thus, it limits the possibility of infection during the treatment. In order to ensure safety, disposable cartridges are packed in blisters and sterilized with ethylene oxide.

When the cartridge is removed from the hand piece, the needle automatically retracts into the adaptor. This prevents from accidental injection and protects against possible infection with viruses that cause serious diseases such as AIDS or viral hepatitis.

The cartridge contains a pigment container. It contains the amount of colour that is enough to perform a standard treatment. The design of the cartridge prevents the pigment from spilling. The shape of the needle blade has been designed in such a way as to obtain the highest precision in placing the colour in the skin, with the least possible harm.



secure.  
CARTRIDGE



made in Poland

[www.bioevolution.com](http://www.bioevolution.com)

# CARTRIDGES ACUPUNCTURE NEEDLES


## secure CARTRIDGE

Innovative cartridges with acupuncture needles and the PROTECT system ensure the highest precision as well as hygiene and safety of the treatment. Extremely sharp flexible needle gently injects the skin without tearing it or causing excessive irritation.

Bleeding is kept to a minimum, so that much more pigment remains in the skin. As a result, the treatment is practically painless. The skin heals very quickly after the treatment.

## RECOMMENDATIONS

- low speed operation at an angle of 90°
- short ejection of the needle
- stretched out skin, delicate injection.

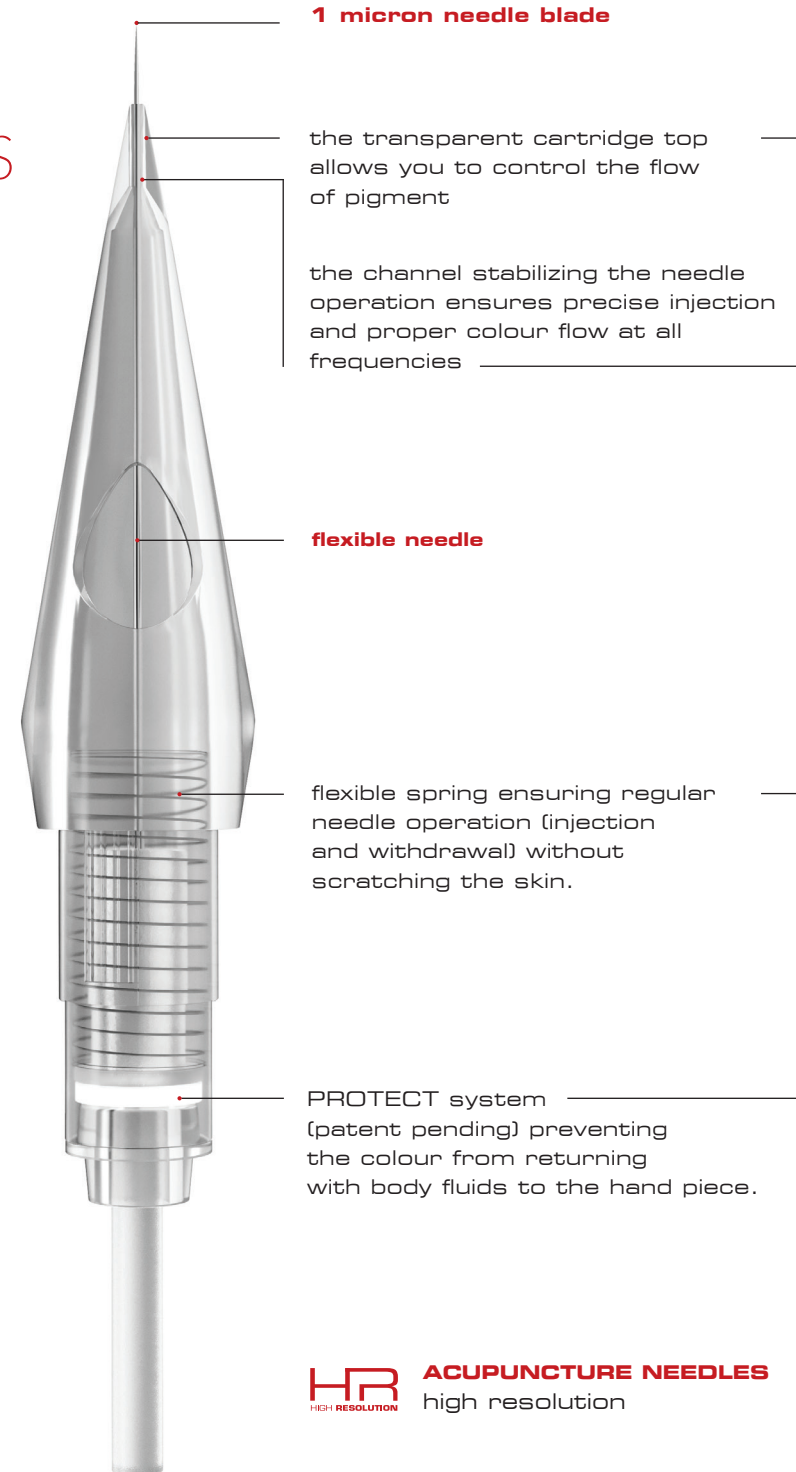

**1 micron needle blade**

the transparent cartridge top allows you to control the flow of pigment

the channel stabilizing the needle operation ensures precise injection and proper colour flow at all frequencies

**flexible needle**

flexible spring ensuring regular needle operation (injection and withdrawal) without scratching the skin.

PROTECT system (patent pending) preventing the colour from returning with body fluids to the hand piece.

3 micron needle blade

rigid needle

## CARTRIDGES TRADITIONAL NEEDLES

secure ●  
CARTRIDGE

Cartridges with traditional needles, thanks to the use of rigid needles are ideal for pigmenting both delicate and hard skins. Properly selected hardness of the needle makes the injection very precise, so we can minimize skin injuries and bleeding. The quick and strong entry of the needle deep into the skin means that we get an even, smooth line of injections. As a result of using such solutions, more pigment is left in the skin and the healing time after the treatment is significantly reduced.

### RECOMMENDATIONS

- operation in any speed range at 90°, 60°, 45°
- short, medium, long needle ejection.

# CARTRIDGES ACUPUNCTURE NEEDLES

secure  
CARTRIDGE

**HR** ACUPUNCTURE NEEDLES  
HIGH RESOLUTION high resolution



1 liner 0.14<sup>HR</sup>



The ultra-precise needle for experienced linergists allows for unprecedented accuracy of highlighting the details and the smallest details of pigmentation.



1 liner 0.16<sup>HR</sup>




The ultra-precise needle for experienced linergists allows for unprecedented accuracy of highlighting details and the smallest details of pigmentation.



1 liner 0.18<sup>HR</sup>



Ultra precise, thin single-point needle for advanced pigmentation techniques, precise lines, hair technique. Designed for experienced linergists.

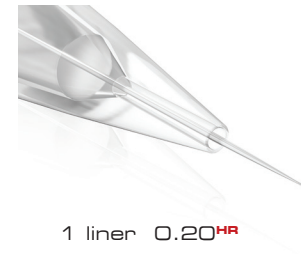

1 liner 0.20<sup>HR</sup>



Very precise needle designed for contours, hair technique, very precise lines, nano technique and scalp pigmentation.



1 liner 0.23<sup>HR</sup>



A single-point needle for precise pigmentation of eyebrows by the imitation of eyebrow hair and powdery method, nano techniques and pixel pigmentation.



# CARTRIDGES ACUPUNCTURE NEEDLES

secure  
CARTRIDGE

**HR** ACUPUNCTURE NEEDLES  
HIGH RESOLUTION high resolution



1 liner 0.25<sup>HR</sup>

A single point needle for pigmentation of precise lines lip contour, nano technique and scalp pigmentation.




1 liner 0.30<sup>HR</sup>

A single-point needle with versatile pigmentation options - it is the most universal needle. Suitable for eyebrow pigmentation, lip contouring, shading as well as pigmentation of eyelid lines and scalp pigmentation. Also suitable for beginner linergists and people who start working with acupuncture needles.



1 liner 0.35<sup>HR</sup>

A single point needle providing a thicker line, for thickening eyelash lines and lines on the eyelid, for pigmenting eyebrows by shading and powdery method. Also recommended for contouring and shading of lips as well as for nipple reconstruction.

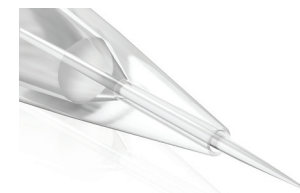

1 liner 0.40<sup>HR</sup>

The needle is intended for clear pigmentation, thick lines, visible contours perfect for eyelid lines, eyebrow shading as well as for nipple reconstruction.

# CARTRIDGES TRADITIONAL NEEDLES

secure  
CARTRIDGE



1 liner 0.20



Very precise needle designed for imitating eye-brow hair, providing very precise lines, nano technique and pigmentation of the scalp and nipple reconstruction.



1 liner 0.25



Single-point needle intended for pigmentation of precise lines, lip contour, nano techniques and hair imitation as well as pigmentation of the scalp and for nipple reconstruction.



1 liner 0.30



Needle designed for making clear lines, contours, eyebrow pigmentation using the powdery and ombre methods. Also recommended for shading and lip contouring, as well as pigmentation of eyelid lines and scalp pigmentation, as well as for nipple reconstruction.



1 liner 0.35




Single-point needle providing a thicker line, for thickening eyelash lines and lines on the eyelids. Also recommended for pigmenting eyebrows using the shading and powdery methods, for contouring and shading the lips, and for the reconstruction of the nipple.



1 liner 0.40



Needle designed for pronounced pigmentation, thick lines, visible contours, perfect for eyelid lines, eyebrow shading as well as for nipple reconstruction.

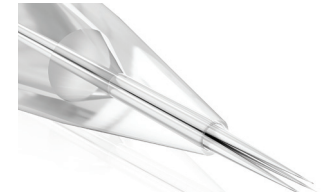


3 liner 0.30



Needle mainly designed for more pronounced eyebrow and lip lines, also suitable for thick skin.

# CARTRIDGES TRADITIONAL NEEDLES


secure  
CARTRIDGE

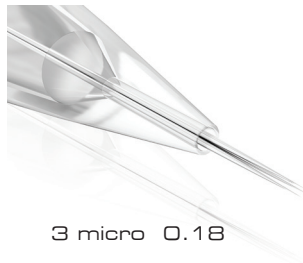

3 micro 0.18




Universal and very popular needle: three precisely connected needles give a beautiful, soft outline of lips shading, powdery eyebrows and eyelid lines.

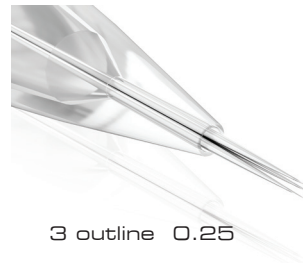

3 outline 0.25




Needle used for pigmentation of eyebrows and lips, for filling and soft shading of larger surface.

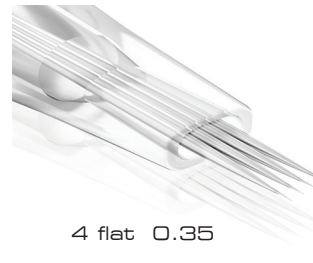

4 flat 0.35




Arranging four needles side by side gives the possibility of faster pigmentation of larger surfaces. The effect of soft powdering and shading of eyebrows and lips, as well as corrective skin pigmentations is obtained .

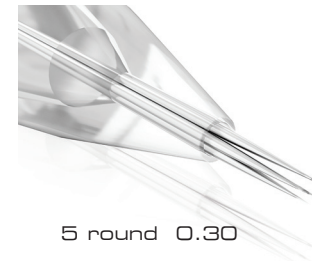

5 round 0.30



It introduces a lot of pigment giving greater colour saturation. Used to fill the lips and shade the eyebrows, it allows you to get a thick line. Not suitable for thin, delicate skin, lip contour or around the eyes.

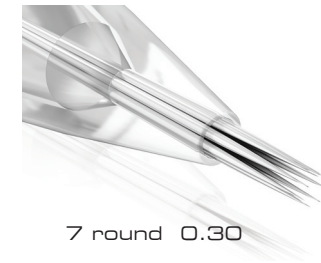

7 round 0.30



It introduces a lot of pigment giving the effect of strong coverage and filling larger areas. Intended mainly for hard skin, resistant to pigment absorbing. Also used for make-up removal.

# EUPHORIA

## COLOURS FOR PERMANENT MAKEUP

BioEvolution EUPHORIA colours for permanent makeup have been developed on the basis of several years of research on the deposition of pigments in the skin. As a result of an in-depth analysis of the composition of colours, types of pigments, anatomy of human skin, it was possible to develop a formula that allows the deposition of more than 90% of pigments contained in the colour, already during the first pigmentation. Thanks to this, the necessity of performing repeated correction treatments has been virtually eliminated.

EUPHORIA colours deposit very well on eyebrows and eyelids. Sensational effects are also achieved with lip pigmentation, which is one of the most difficult areas for pigmentation. The lips are very well saturated with pigment red colour already after the first treatment, which perfectly covers the natural colour of the lips.

By pigmenting with EUPHORIA colours you can enjoy expressive, vivid colours and a beautiful, natural effect. EUPHORIA colours are produced in accordance with the Cosmetics Directive of the European Union (EU resolution ResAP 2008/1).

capacity: 10 ml

# COLOUR CHART

## EYEBROWS

BROWN BLACK 210	LAVA 230	TOFFE 250
CHOCOLATE 220	DARK BROWN 235	KHAKI 255
MOCCA 215	HAZELNUT 240	BLONDE 260
GRAPHITE 225	GOLDEN BROWN 245	

## LIPS

LIGHT ROSE 310	CORAL 330	NUDE 350
RASPBERRY 315	TEA ROSE 335	PEACH 355
RED WINE 320	RED BROWN 340	CARAMEL 360
CHERRY 325	TERRACCOTTA 345	SALMON 365

## EYELID LINES

BLACK 100	MEDIUM GREEN 061
DEEP BLUE 050	DARK GREEN 062
AZURE BLUE 051	
LIGHT GREEN 060	

## CAMOUFLAGE

WHITE 010
CAMO 020
YELLOW 030
ORANGE 040



## COLOURS FOR PERMANENT MAKEUP

The new line of BioEvolution QLINE colours is based on the latest technology, which limits the pigmentation of eyebrows, eyes and lips to unwanted colours to the maximum.

The BioEvolution QLINE colour formula has been developed to make the make-up look natural. The make-up colours are soft, delicate, with a subtle satin effect. At the same time, the colours are very saturated and vivid.

Thanks to the unique technology of small molecule pigments, BioEvolution QLINE colours have unique covering properties. They settle very easily and quickly in the skin.

Produced in accordance with the Cosmetics Directive of the European Union (EU resolution ResAP 2008/1).



capacity: 5 ml

# COLOUR CHART

## EYEBROWS

BROWN 800	BROWN 810	BROWN 820 MX	BROWN 840
BROWN 800 MX	BROWN 815	BROWN 825	
BROWN 805	BROWN 815 MX	BROWN 830	
BROWN 805 MX	BROWN 820	BROWN 835	

## LIPS

RED 700	RED 720	RED 740	RED 760
RED 705	RED 725	RED 745	
RED 710	RED 730N	RED 750	
RED 715	RED 735	RED 755	

## EYELID LINES

BLACK 810 EL	BLUE 110
BLACK 120	BLUE 120
GREEN 110	BLUE 130
GREEN 120	

## MODIFIERS

ORANGE 001
WARM UP
RED STOP
GOLDEN BROWN

\* MX - Colour with a higher content of orange pigment - reduces the risk of „gray eyebrows” effect.

\*\* EL - Pigment only for eyelid lines.

### IMPORTANT NOTICE!

Do not mix QLINE pigments with pigments not belonging to this line.

# MICROBLADING

The BioEvolution Qline Microblading colours are intended for eyebrow pigmentation using the microblading method and other manual techniques. Colours are characterized by a creamy consistency with an ideal density. They settle in the skin easily and provide a long-lasting effect without discoloration.



capacity: 10 ml

## BRWI

DIAMOND BROWN	MEDIUM BROWN	COLD SKIN BROWN	DARK BROWN	SUN KISSED	WARM SKIN BROWN
					
BLACK AMBER	BLACK ONYX	SKIN	ORANGE	RED STOP	
					

Notice:

Depending on the microblading technique used and the depth at which the colour is introduced, you can get a cooler shade of hair after healing.